


# NEOFUSER ONCO

## MICROINFUSION PUMPS WITH SINGLE FLOW RATE FOR 5-FU

A continuous portable infusion pump, specifically designed for oncology needs. Pumps with preset flow rate, suitable for application in regional anesthesia, wound infiltration, as well as in the post-operative pain therapy, among others.

- Pump volume: 100(M), 275(L), 450(X) ml.
- Basal rates: 0.5, 1.0, 1.5, 2.0, 2.5, 3.0, 4.0, 5.0, 8, 10, 20, 50, 100, 200, 250 ml/h.
- Protection against UV light up to 520 nm.
- Transparent material allows for detection of air bubbles and particles.
- Not made with latex, BPA, or DEHP.
- Luer-Lock patient connector with color coded flow rate control unit.
- Efficient particle and self-venting 0,2 µm filter.
- Pressure-Activated Valve for one-way drug delivery.
- The clamp is used to start and stop the flow.

Basal Flow Rate, ml/hr	Color code	Ref. Nominal Volume		
		100 ml	275 ml	450 ml
0,5	■	M05-O	L05-O	-
1	■	M10-O	-	-
1.5	■	-	L15-O	-
2	■	M20-O	L20-O	X20-O
2.5	■	M25-O	-	-
3.0	■	-	-	X30-O
4	■	M40-O	L40-O	X40-O
5	■	M50-O	L50-O	X50-O
8	■	-	-	X80-O
10	■	-	L100-O	X100-O
20	□	-	-	X200-O
50	■	-	-	X500-O
100	■	-	-	X1000-O
200	■	-	-	X2000-O
250	■	-	-	X2500-O

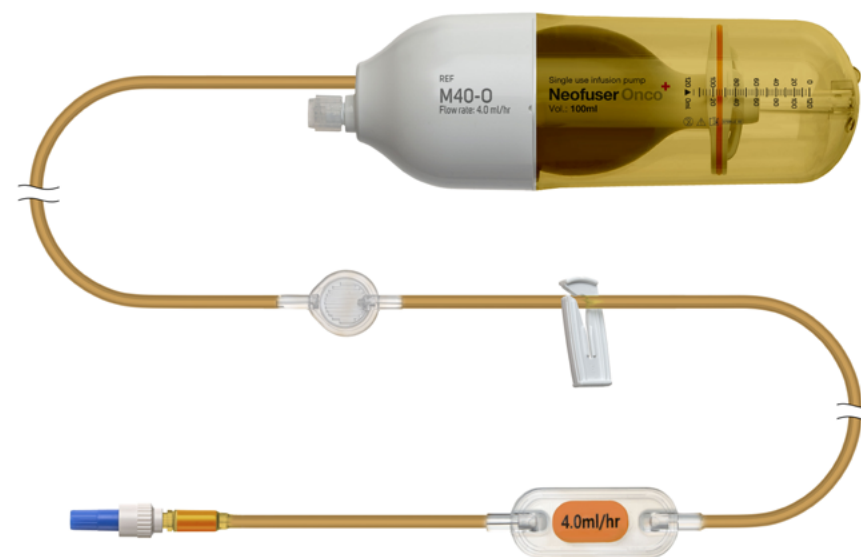

Find out more at [www.invumed.com](http://www.invumed.com)

INVUMED Innovations LTD. Cyprus. 23 John Kennedy Avenue «Globe House» - 1075 Nicosia  
invumed.com / info@invumed.com / +357 26220591  
© 2022 INVUMED Innovations LTD. All rights reserved. 22.12.22